


Ambu® VivaSight™ 2 SLT

The second generation of Ambu endotracheal tubes for One Lung Ventilation

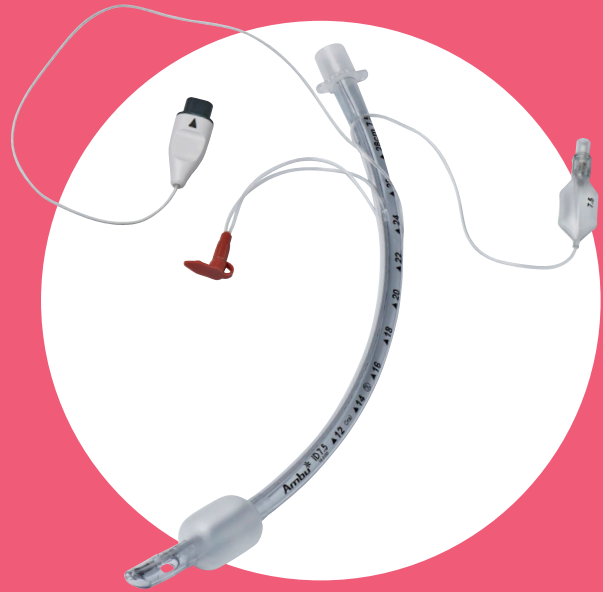

KEY BENEFITS

Facilitates endotracheal intubation

Provides continuous visualization of tube and endobronchial blocker position

Significantly reduces the need to use a bronchoscope

Compatible with Ambu® aView™ 2 Advance or Ambu® aView™ Monitor

DATASHEET

VivaSight 2 SLT is a sterile, single use endotracheal tube with integrated camera and is intended for:

- Oral intubation procedures.
- General inspection of the airways and for visualization during intubation procedures.
- Temporary artificial airway in adults requiring mechanical ventilation.

VivaSight 2 SLT is similar in design to conventional Single Lumen Tubes, but has a built-in video camera that connects to an Ambu displaying unit, making it possible to monitor tube placement in real-time.

Real-time monitoring during One Lung Ventilation

Due to the integrated camera and continuous visualization, VivaSight 2 SLT:

- Facilitates endotracheal intubation.
- Provides continuous visualization of endobronchial blocker position.
- Ideal for monitoring of correct tube placement.
- Significantly reduces the need to use a bronchoscope.

The full Ambu OLV portfolio includes VivaSight 2 DLT, VivaSight™ 2 SLT, Ambu® VivaSight™ Endobronchial Blocker Tube, Ambu® aScope™ 4 Broncho Slim and other Ambu bronchoscopes with a small diameter. These bronchoscopes, along with VivaSight 2 DLT and VivaSight 2 SLT, are all compatible our latest generation portable full-HD endoscopy system, Ambu® aView™ 2 Advance.

VivaSight 2 SLT can be connected with Ambu's full HD monitor, aView 2 Advance, for high quality imaging and Dual View.

Ambu

SPECIFICATIONS

Ambu® VivaSight™ 2 SLT	SL 7.0 mm	SL 7.5 mm	SL 8.0 mm
Size	7.0 mm	7.5 mm	8.0 mm
ID - Effective Inside Diameter (mm)	4.6	4.8	5.2
OD (mm)	10.0	10.5	11.0
Length (±3 mm)	323	323	323

General

Usage: Single-use only

Pack: Supplied individually packaged, ETO-sterilized

Shelf life: 36 months

Storage: Store at temperatures between 10 °C - 25 °C in a dry, cool and dark place

Materials

Tube: Medical-grade PVC; additional injection port for rinsing camera optics

Cuff: Medical-grade PVC, high volume, low pressure

Tube hardness: 85 (Shore A)

This product is not made with phthalates

Camera

Image sensor: CMOS

Light Source

LEDs: 2 white

ORDERING SPECIFICATIONS

Ambu® VivaSight™ 2 SLT Kits	Order number	Quantity/box
Ambu® VivaSight™ 2 SLT - 7.0 mm (with adapter cable)	412701000	5
Ambu® VivaSight™ 2 SLT - 7.5 mm (with adapter cable)	412751000	5
Ambu® VivaSight™ 2 SLT - 8.0 mm (with adapter cable)	412801000	5

Ambu® VivaSight™ 2 SLT + EBB Kits	Order number	Quantity/box
Ambu® VivaSight™ 2 SLT + EBB Kit - 7.0 mm, 9 Fr (with adapter cable)	412709000	5
Ambu® VivaSight™ 2 SLT + EBB Kit - 7.5 mm, 9 Fr (with adapter cable)	412759000	5
Ambu® VivaSight™ 2 SLT + EBB Kit - 8.0 mm, 9 Fr (with adapter cable)	412809000	5

Monitors

Ambu® aView™ 2 Advance Monitor*	405011000	1
---------------------------------	-----------	---

* Monitors item number is different in each country. Please contact your local sales rep for the right number in your location.



Ambu® VivaSight™ 2 SLT



Ambu® VivaSight™ 2 SLT connected to Ambu® aView™ 2 Advance Monitor

Ambu

Ambu A/S
 Baltorpbakken 13
 2750 Ballerup
 Denmark
 T +45 72 25 20 00
 F +45 72 25 20 50
 ambu.com

Rx ONLY

